

PICC access with power injectability.


The SYNERGY™ CT PICC offers the ease and effectiveness of PICC access with power injectability, creating a SYNERGY of intravenous therapy options for healthcare practitioners.

Polyurethane catheter construction provides radiopacity, durability and strength

Low profile hub allows for suture or sutureless securement

Easy legible identification on clamps and extensions, indicating power injectability and maximum pressure

Simple grip luers are easy on gloves and maintain strength and crack resistance

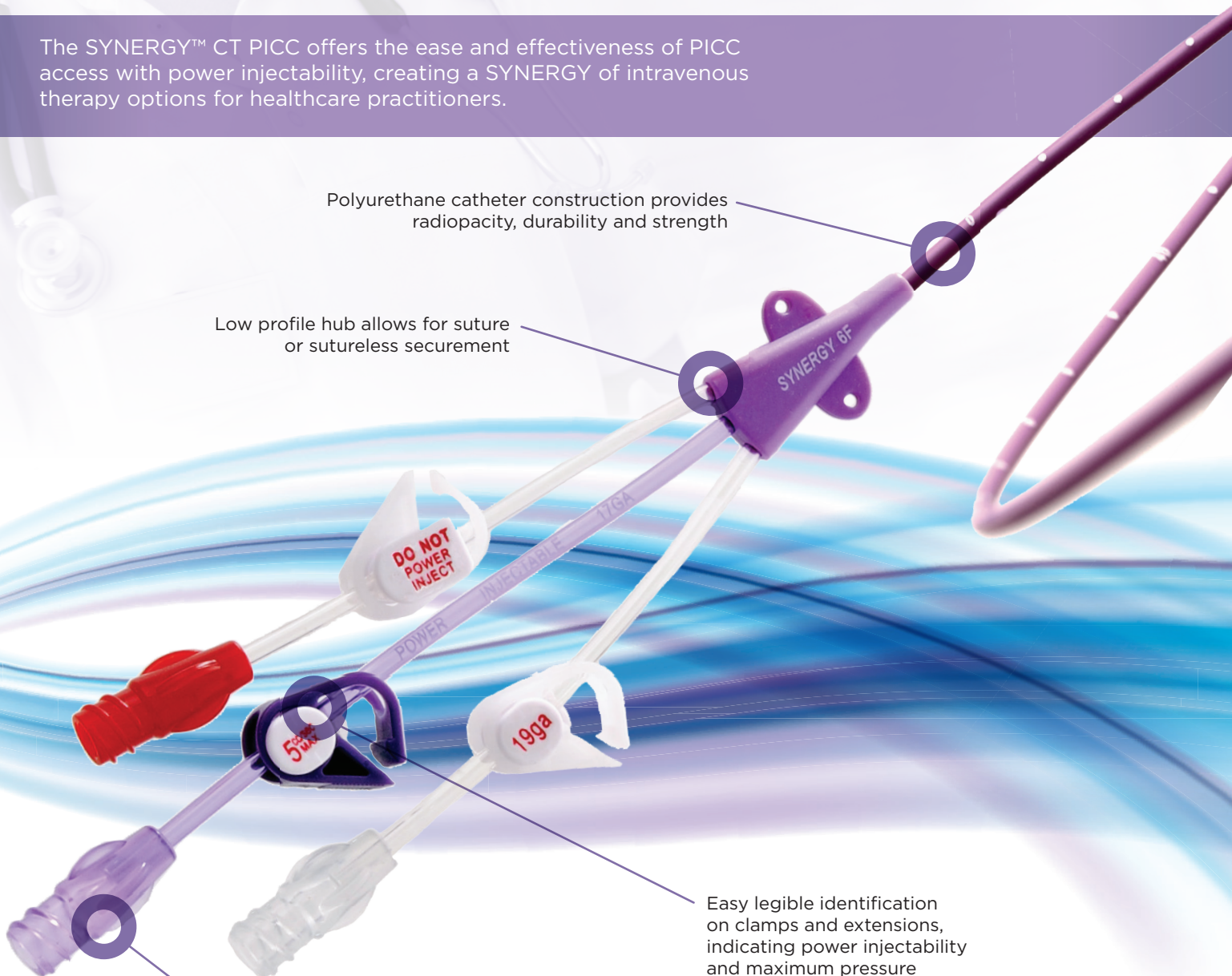

Features and Benefits

Maximum Injection Rates – allow for the injection of contrast media at rates up to 5ml/sec for CECT scans.

CVP Monitoring – SYNERGY power PICCs are indicated for central vein pressure monitoring.

Polyurethane Construction – reverse taper technology helps reduce post insertion bleeding while maintaining radiopacity and patient comfort.

Standardized Trays and Components – Trays are organized with familiar industry components for clinician safety and user friendly procedural steps ; including tearaway sheaths, echogenic safety needles, and premium guidewires.

Multiple Kit Configurations – SYNERGY kits are available in Nursing, Basic, I.R. Long Wire, and Max Barrier tray configurations.

Low Profile Hub – allows for better patient comfort and easy securement.

SYNERGY Size Configurations

| Single Lumen | Double Lumen | Triple Lumen |
|--------------|--------------|--------------|
| 4fr | 5fr | 6fr |
| 5fr | 6fr | |

Standard Tray Components

| | Quantity |
|---|------------|
| SYNERGY™ CT PICC | 1 |
| Needle Free Valve | 1, 2, or 3 |
| 0.15 x 75cm Twisted Wire Stylet | 1 |
| .030 I.D. Adaptor W/Sideport | 1 |
| 21GA X 7CM SAFETY SELD. ECHO.NEEDLE | 1 |
| PTFE Tearaway Introducer | 1 |
| Nitinol Straight Mandrel GW W/ Tungsten Tip | 1 |
| 10cc Disposable Luer Lock Syr | 1 |
| Safety Scalpel #11 | 1 |
| 24" Paper Tape Measure | 1 |
| PICC Stylet TAG | 1 |
| PICC Patient Info Guide | 1 |
| PICC Chart Sticker | 1 |
| PICC IFU | 1 |
| Securement Device | 1 |

For ordering information or for question about the Health Line products, please contact your Health Line representative or send us an email at info@hlic.net.



Health Line International Corp.
803 North 1250 West #1
Centerville, Utah 84014 USA
Phone: (801) 773-7798
Fax: (801) 820-8007
www.hlic.net

ISO13485
Certified Quality System